


# C18 UH2

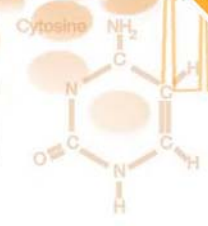

## Filter & Guard Options

- Guard system for all 3µm, 5µm and 10µm phases
- Low volume in-line filters for LC and UHPLC
- Maintain chromatographic integrity

Fortis Guards and filters are designed to ensure that erroneous materials do not find their way onto the more important and expensive analytical column. Guards are available in sizes to match all analytical and preparative column dimensions. Filters are particularly suitable for short fast LC/MS (Pace™) columns and UHPLC columns.



- Guard system for all 3µm, 5µm and 10µm phases
- Quick replacement cartridges
- Highly Cost Effective



- Preparative guard system 10mm & 21.2mm
- Quick replacement cartridges
- Highly Cost Effective



- Filter for all LC columns
- Low volume In-line filters
- Change over time - seconds not minutes

Fortis in-line filters are fingertight direct connect design, fitting in between the column and the conventional peek fitting to filter out particulate matter, it contains low dead volume and pressure. In-line filters are ideal for very short fast columns such as Fortis Pace™ LC/MS columns where extra packed bed from a guard would be detrimental.

In-line filters are also available in UHPLC format, capable of withstanding the elevated pressures involved.

3µm Fortis C18 Cartridges	
F18-040003G/2	10x4mm Fortis C18 3µm Guard pk 2
F18-040003G/4	10x4mm Fortis C18 3µm Guard pk 4
F18-020003G/2	10x2mm Fortis C18 3µm Guard pk 2
F18-020003G/4	10x2mm Fortis C18 3µm Guard pk 4

5µm Fortis C18 Cartridges	
F18-040005G/2	10x4mm Fortis C18 5µm Guard pk 2
F18-040005G/4	10x4mm Fortis C18 5µm Guard pk 4
F18-020005G/2	10x2mm Fortis C18 5µm Guard pk 2
F18-020005G/4	10x2mm Fortis C18 5µm Guard pk 4

In-line Filters	
2-SAV5	2µm In-line filter pk 5
2-SAV10	2µm In-line filter pk 10
UHPSAV2	UHPLC In-line filter pk 2
UHPSAV4	UHPLC In-line filter pk 4